



Vscan Air™

Ultrasound



Authorized Distributor
GE Healthcare

Vscan Air™

See More. Treat Faster.

Vscan Air handheld ultrasound probe is the 4th generation in the Vscan family. A truly wireless, lightweight handheld ultrasound system that enables the flexible performance of whole-patient scanning anywhere from the convenience of your pocket.

Probe Features

- Dual-probe technology
- Ergonomic design
- Pairing across multiple devices
- IP67 waterproof rating and compatibility with high-level disinfection techniques
- MIL-810G military-drop standards
- Up to 50-minute scan time
- Recharge within 75 minutes (from 10% to 90%)
- Use any Qi-compliant wireless charger


Connection App

- Intuitive user interface
- Minimized number of keys and intuitive thumb-controllable touchscreen user interface
- Single key/gesture to control freeze/unfreeze store, color on/off, and gain/depth control
- Configuration of probe button (freeze or store)
- Easily accessible presets — available without scrolling
- User-selectable default preset


Small and Lightweight





Apple iPhone and iPad

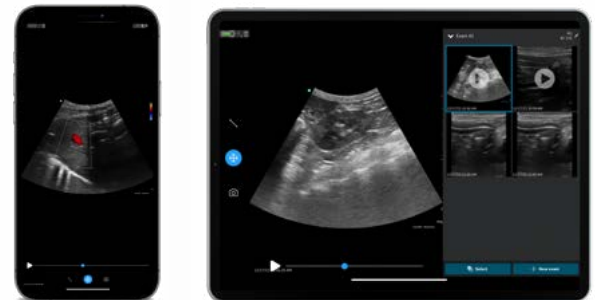

Android Smart Phone and Tablet

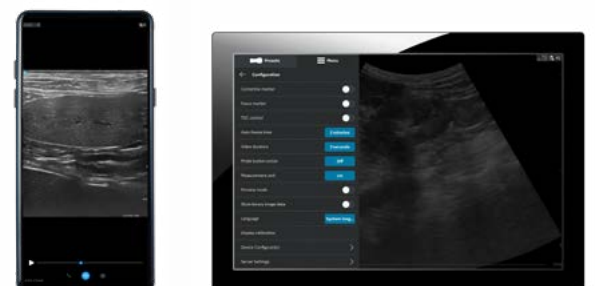

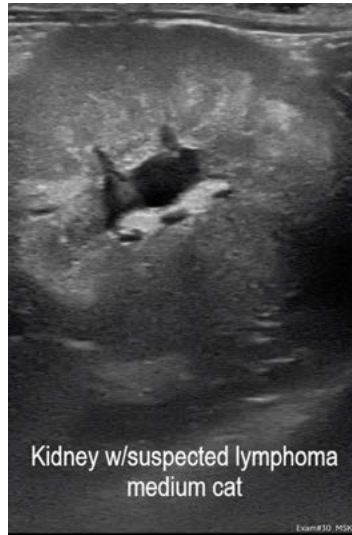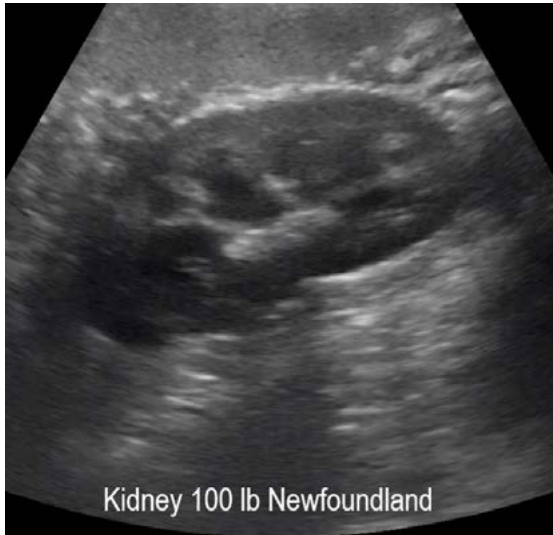

SignalMAX™ Technology

Historically, high-quality images have required large and expensive hardware. Additionally, handheld ultrasound devices, though small and ergonomic, were lower image quality.

With SignalMAX technology, you don't have to compromise on image quality to have a probe that comfortably fits in your hand.

Vscan Air's dual-probe uses the same piezoelectric transducer technology as GE Healthcare's leading ultrasound brands. Each side of the probe is tailor-made for the requested task.

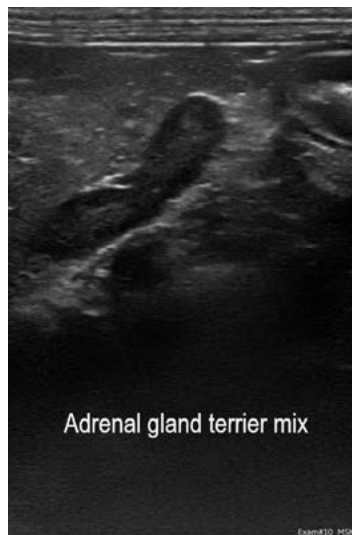


Image Storage & Access

DICOM®, including query from modality wordlist server, store to DICOM PACS, receive storage commitment



JPG/MPG export to shared network folders (supports importing into EMR)

Share auto-anonymized images with other apps

No cloud required — easily share images with owners, colleagues, and social media

Secure DICOM including the option of peer authentication

For Here. There. Anywhere.

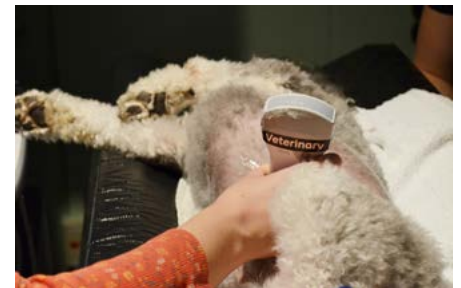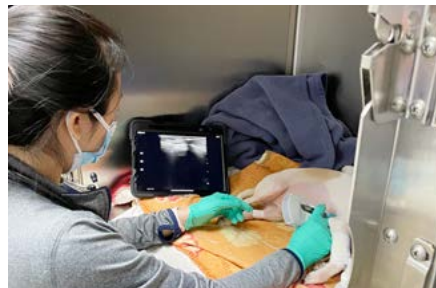

**To learn more and how
Antech can benefit your
veterinary organization**

VISIT
antechdiagnostics.com

CALL
1-800-268-5354



ANTECH

Sound Imaging 

Antech United States

1-800-872-1001

Antech Canada

1-800-341-3440

Antech Sound Imaging

1-800-268-5354

Antech Sound Imaging Support

1-800-819-5538

Antech Sound Equine

1-800-955-6353

antechdiagnostics.com



© 2022 Antech Diagnostics, Inc. All rights reserved. Antech, the Antech logo, and all other trademarks used herein are the registered trademarks of Antech Diagnostics, Inc. or its affiliates. 042122